

Maxso[®]

X92


Implant handpiece

Experience unrivaled cutting power with the X92 implant handpiece by Beyes. This powerful unit features both internal and external irrigation for maximum cooling where it is needed the most. Its outstanding ergonomic grip offers the perfect balance, resulting in high precision operation. A titanium-coated, scratch-resistant surface ensures longstanding durability. Beyes' X92 implant handpiece has been engineered to ease the most challenging procedures.



Push Button Chuck

Effectively cools the entire surface of the bur that avoids overheating and provides better patient comfort.

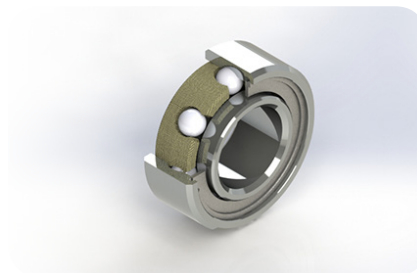


Ceramic Ball Bearings

Sturdier yet lighter than steel, ceramic bearings increase both efficiency and longevity of the unit.



Dual Irrigation

An internal and external irrigation system supplies ample cooling water to the operating area.

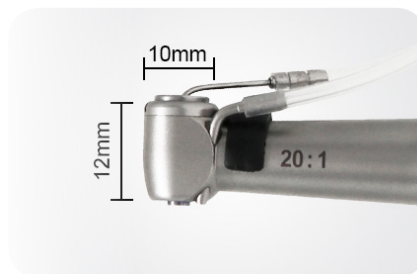

Mini-sized Head

Experience unobstructed access with an extra small head size.

Key Specification

Gear Ratio	20:1 Reduction
Dual Irrigation	External
Body Material	Stainless Steel
Max Speed	2,000 min-1
Max Torque	80 Ncm
Sterilization	Autoclavable up to 135°C

Order Information



ISO E-Type Connection

Code: SL3019

Model: X92

Description: Contra Angle, 20:1 Implantology, Non-Optic, Push Button, CA Burs



Beyes[®]

The Next Standard

Beyes Dental Canada Inc.
23-595 Middlefield Road
Toronto, Ontario, M1V 3S2
Canada

Tel: 1-416-281-1888

Email: info@beyes.ca

Website: www.beyes.ca

Your Dealer for Beyes Products

Scan the QR code for more information on X92



Printed in Canada
Document #: ENB019
Rev.1 /12.07.19

© Beyes Canada Inc. 2019. All Rights Reserved.
Specifications are subject to change without notice.
All brands are holders of their respective trademarks.
All Beyes products are manufactured in strict compliance with Health Canada regulations ensuring quality.
Manufacturing facility is ISO 13485 certified under Canadian Medical Devices Conformity Assessment System (CMDCAS).