


UP TO TWO MM

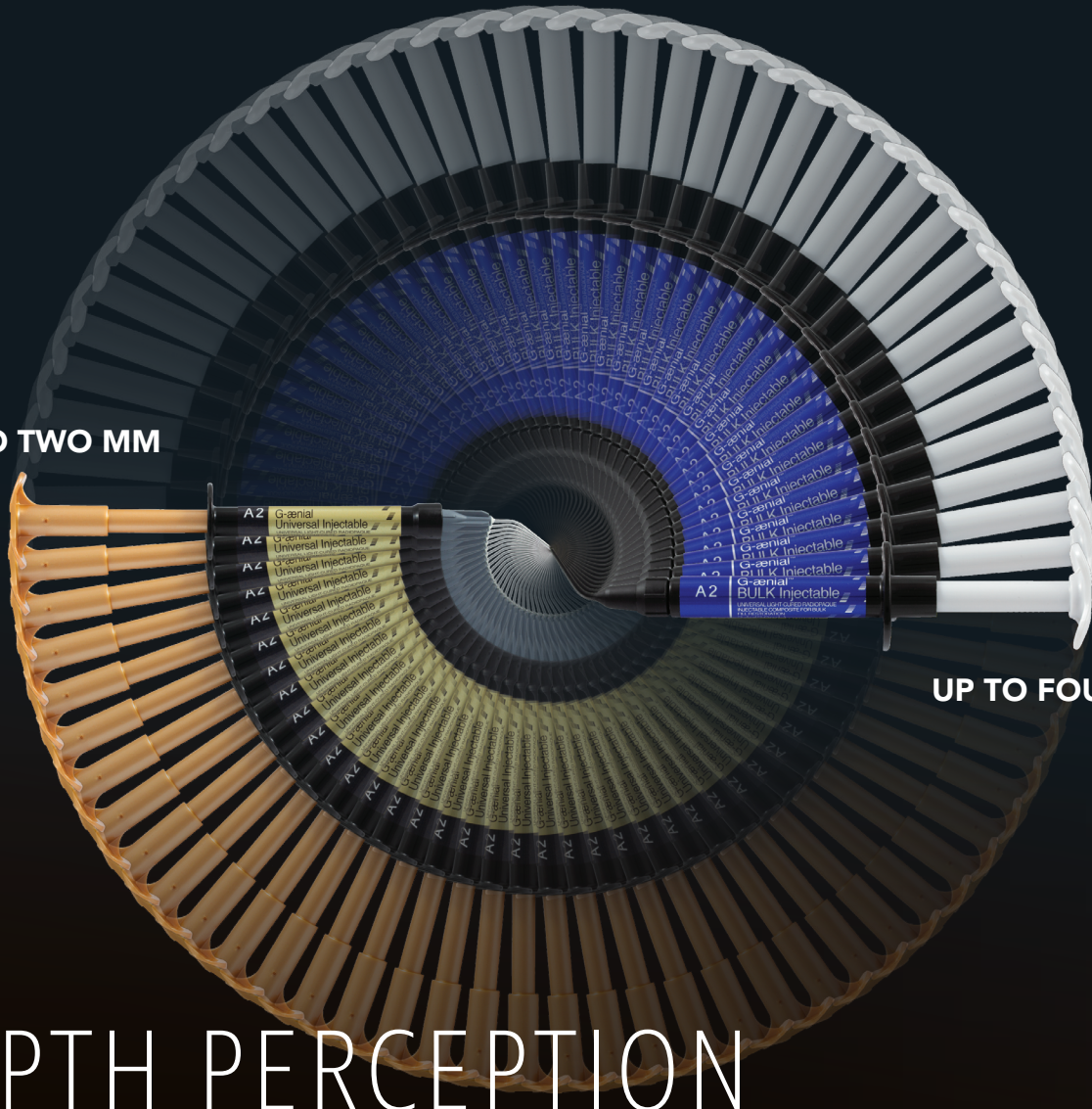

UP TO FOUR MM

DEPTH PERCEPTION

G-ænial™ BULK Injectable
G-ænial™ Universal Injectable
G-ænial™ Flo X
Family of Injectable Composite Restoratives





REDEFINE THE WAY YOU WORK WITH GC INJECTABLES!

G-ænial™ composites are great for all cavity depths and classes!

Class	Cavity Depth	Product Recommendation	
I / II	Up to 2 mm	G-ænial™ Flo X or	
		G-ænial™ Universal Injectable	
I / II	Up to 4 mm	G-ænial™ BULK Injectable	
I / II	4-6 mm	G-ænial™ BULK Injectable & G-ænial Sculpt®	
I / II	6 mm +	G-ænial™ BULK Injectable & G-ænial Sculpt®	

When used with G-ænial Sculpt® and G-Premio BOND™, you can offer your patients the ultimate combination of strength and esthetics.





FOR SMALLER POSTERIOR RESTORATIONS...

G-ænial™ Flo X



Light-Cured Radiopaque Flowable Composite

- Homogenously dispersed ultra-fine barium-glass particles provide excellent handling
- Optimal viscosity for lining; flows easily and ensures excellent adaptation to cavity walls
- High radiopacity, easily detectable on radiograph even at low thickness




Featuring high radiopacity, easily detectable on radiograph even at low thickness, excellent handling properties, and optimal durability, G-ænial™ Flo X is the perfect partner to G-ænial™ Universal Injectable.

- Indicated for use in small restorations and as a cavity liner
- Highly radiopaque and highly flowable
- Provides easy placement and flow into preparations for excellent handling

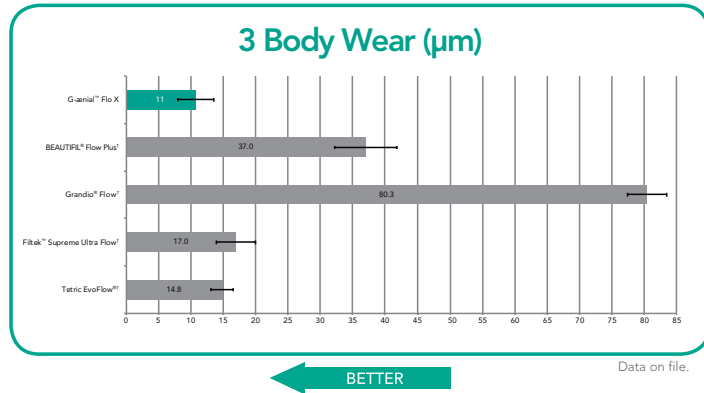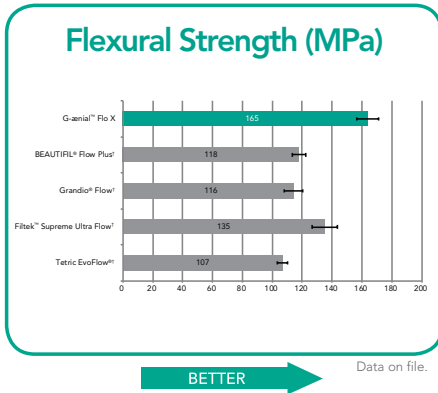

Lining or Base

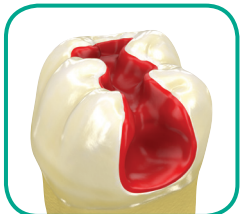


Fissure Sealing



Small Marginal Defect



Filling Undercuts

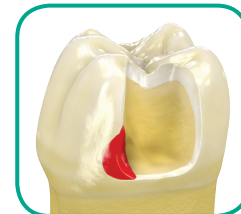

Tunnel Restorations



FOR EXCEPTIONALLY STRONG RESTORATIONS...

G-aenial™ Universal Injectable

Injectable High-Strength Ultra-Fine Particle Composite

- Injectable, shapeable paste
- Virtually no slumping, stays where you put it
- Create anatomy as you inject
- Excellent adaptation to cavity preparation
- Fast and easy placement, self polishing in the mouth
- High radiopacity, easy to identify on x-ray




An injectable high-strength ultra-fine particle composite with ideal viscosity, handling and adaptation characteristics that may be used for long-lasting esthetic restorations.

The homogeneously dispersed ultra-fine barium-glass particles formula (150nm particle size, 69% wt filler loading) provides high flexural strength and wear resistance for a durable and long lasting restoration.






Metal tip may be bent with fingers if necessary to better access posterior distal region.

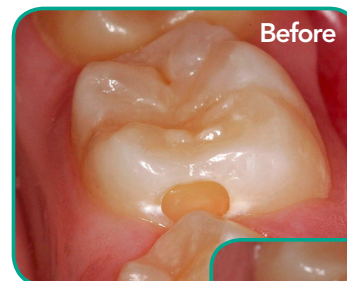


Courtesy of Robert Marus, DDS, Yardley, PA

G-aenial™ Universal Injectable features high flexural strength, for protection against chipping and fracture, and outstanding wear resistance for excellent durability and gloss retention.

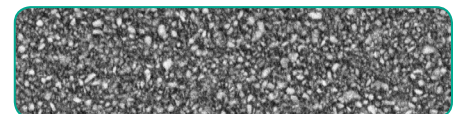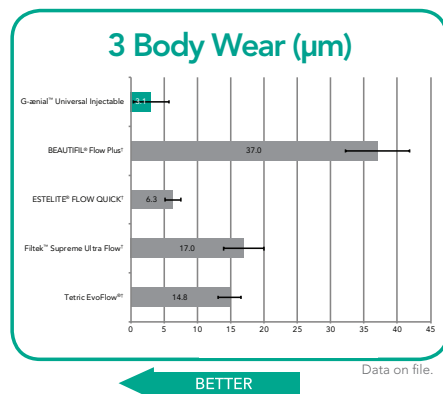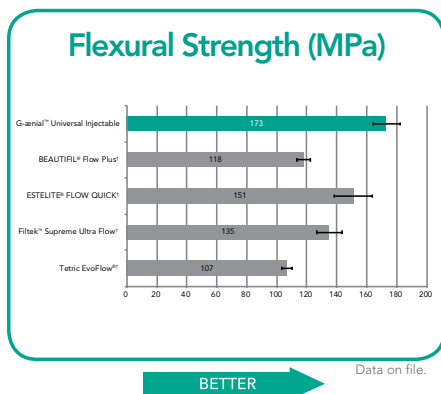

G-aenial™ Universal Injectable
Ultra-Fine Particles

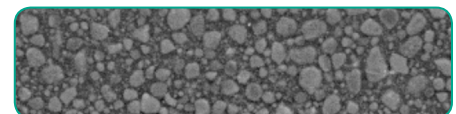

Competitive Product
Larger Particles

† BEAUTIFIL® Flow Plus is a registered trademark of Shofu Dental Corporation. ESTELITE® Flow Quick is a registered trademark of Tokuyama Dental Corporation. Filtek™ Supreme Ultra Flow is a trademark of 3M. Tetric EvoFlow® is a registered trademark of Ivoclar Vivadent.



FOR FAST, EFFICIENT BULK FILL RESTORATIONS...


G-ænial™ BULK Injectable Injectable High-Strength Ultra-Fine Particle Composite

- The ideal convergence of adaptation, strength and esthetics
- Highly shapeable; stays where you put it
- Fill up to 4mm, no capping or veneering required
- Excellent adaptation to cavity preparation



A restorative composite that is ideal for bulk filling to the occlusal surface of large cavities, by virtue of its injectable viscosity and excellent physical properties even superior to paste type composites*.



- Injectable composite with excellent viscosity and thixotropy for ideal bulk filling
- Can achieve highly durable restoration via high wear resistance and flexural strength (unlike competitive composite*)
- Excellent polishability and gloss retention with self-polishing ability
- Guaranteed depth of cure of 4mm which allows for bulk filling in one layer

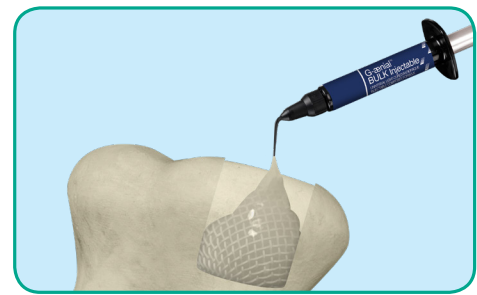

**The Ideal Convergence of
Adaptation, Strength & Esthetics**

G-ænial™ BULK Injectable provides the right combination of strength and esthetics to provide a durable restoration that looks great at every recall!

- High flexural strength to protect chipping or fracture
- Outstanding wear resistance to provide ultimate durability, with excellent color and gloss retention
- Low volumetric shrinkage and low shrinkage stress to preserve marginal integrity and protect against recurrent caries

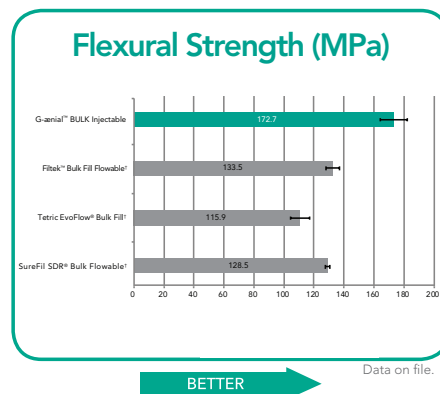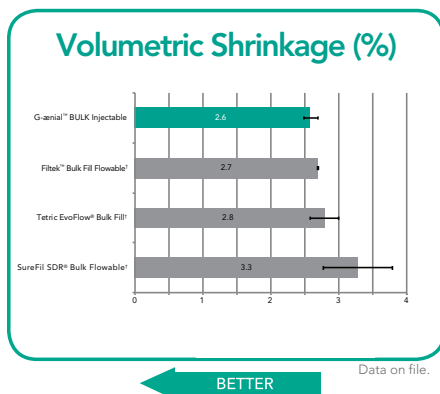

† Filtek™ Bulk Fill Flowable is a trademark of 3M. Tetric EvoFlow® Bulk Fill is a registered trademark of Ivoclar Vivadent. SureFil SDR® is a registered trademark of Dentsply Sirona. *Data on file.



FAST, EASY TO USE & PROVIDES GREAT RESULTS!



G-aenial™ Flo X

G-aenial™ Flo X Refills

Contains: One 3.6g (2.0mL) syringe, 20 Dispensing Tips III (plastic) and one light protective cap.

008317 A1	008318 A2	008319 A3
008320 A3.5	008321 A4	008322 AO2
008323 AO3	008324 CV	

G-aenial™ Universal Injectable

G-aenial™ Universal Injectable Refills

Contains: Two 1.7g (1.0mL) syringes, 10 long dispensing tips and 10 short dispensing tips.

012364 A1	012365 A2	012366 A3
012367 A3.5	012368 A4	012369 B1
012370 B2	012371 CV	012372 CVD
012373 AO1	012374 AO2	012375 AO3
012376 JE	012377 AE	012378 XBW
012379 BW		

G-aenial™ BULK Injectable

G-aenial™ BULK Injectable Refills

Contains: One 3.4g (2.0mL) syringe, 20 Dispensing Tips III (plastic) and one light protective cap.

012394 A1	012395 A2
-----------	-----------

Complementary Products

G-aenial Sculpt®

G-aenial Sculpt® Syringe Refills

Contains: One 2.0mL Syringe.

009152 A1	009153 A2	009154 A3
009155 A3.5	009156 A4	009168 AE
009161 AO1	009162 AO2	009163 AO3
009157 B1	009158 B2	009159 B3
009164 BW	009160 C3	009166 CV
009167 JE	009165 XBW	

G-aenial Sculpt® 20-Unitip Refills

Contains: 20 0.16 mL unitips.

009169 A1	009170 A2	009171 A3
009172 A3.5	009173 A4	009174 B1
009175 B2	009176 B3	009177 C3

G-aenial Sculpt® 10-Unitip Refills

Contains: 10 0.16mL unitips.

009185 AE	009178 AO1	009179 AO2
009180 AO3	009180 BW	009183 CV
009184 JE	009182 XBW	

G-Premio BOND™

009275 G-Premio BOND™ Kit

Contains: One 5mL bottle, 50 disposable applicators (fine), 20 disposable dispensing dishes.

009277 G-Premio BOND™ 5mL Bottle Refill 009276 G-Premio BOND™ Unit Dose Package

Contains: 50 unit doses (0.1 mL each).