


# APEXPILOT® G2


## Cordless Endo Motor

ApexPilot G2 can be connected to the matched apex locator to realize apex locating function. By connecting the ApexPilot G2 to the matched apex locator, the position of the file tip inside the canal can be monitored during the procedure and many automatic functions such as auto apical reverse can be activated.

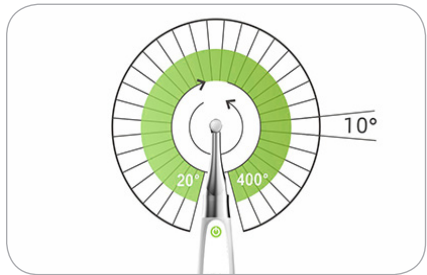

### High Performance Brushless Motor

Adjustable reciprocating angle (10° interval; adjustment range: 20°-400°); Accurate rotating angle reduces risk of instrument fracture and enables higher preparing efficiency.

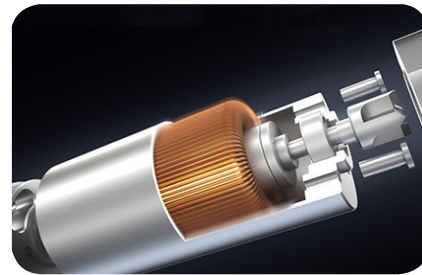

### Efficient & Safe

Super high torque, high speed, fast response, more accurate output, more stable, to ensure more efficient & safe preparation.



### Lower Noise

Smooth rotation, lower noise and small vibration in handpiece. Doctors feel good during operations. Patients feel more comfortable during operation.



### 6:1 Super Mini Contra-angle

This 360° rotatable and super mini contra-angle allows clearer operation field.

## Key Specifications

<b>Speed Range</b>	100-1200RPM
<b>Torque Range</b>	0.4-5.0N.cm
<b>Lithium Battery</b>	3.7V /2000mAh
<b>Power Input</b>	~100V-240V 50Hz/60Hz 0.5-0.2A
<b>Power Output</b>	DC5V/1A
<b>Wireless charging Frequency Range</b>	112-205KHz
<b>Environment Temperature</b>	+5°C~ +40°C
<b>Relative Humidity</b>	30% ~ 75%
<b>Atmospheric Pressure</b>	70kPa ~ 106kPa
<b>Warranty</b>	12 Months

## Order Information

**Code:** AL2060

**Model:** ApexPilot G2



**Beyes Dental Canada Inc.**  
23-595 Middlefield Road  
Toronto, Ontario, M1B 3M5  
Canada  
**Tel:** 1-855-603-1888  
**Email:** info@beyes.ca  
**Web:** www.beyes.ca



**Lotus NL B.V.**  
Koningin Julianaplein 10, le Verd,  
2595AA, The Hague,  
Netherlands  
**Tel:** +31645171879 (English)  
+31626669008 (Dutch)  
**Email:** info@beyes.ca

### Your Dealer for Beyes Products

Scan the QR code  
for more products



Printed in Canada  
File #: ENB031  
Rev.1/04.02.21

© Beyes Canada Inc 2021. All Rights Reserved.  
Specifications are subject to change without notice.  
All brands are holders of their respective trademarks.  
All Beyes products are manufactured in strict compliance  
with Health Canada regulations ensuring quality.  
Manufacturing facility is ISO 13485 certified under Canadian  
Medical Devices Conformity Assessment System (CMDCAS).