




MOJAVE

# MOJAVE<sup>®</sup>

Superior Performance.  
Energy Efficiency.  
All In A Singular Dry Vacuum.

**AIR**  
**TECHNIQUES** equipped for life<sup>®</sup>

# Superior Evacuation

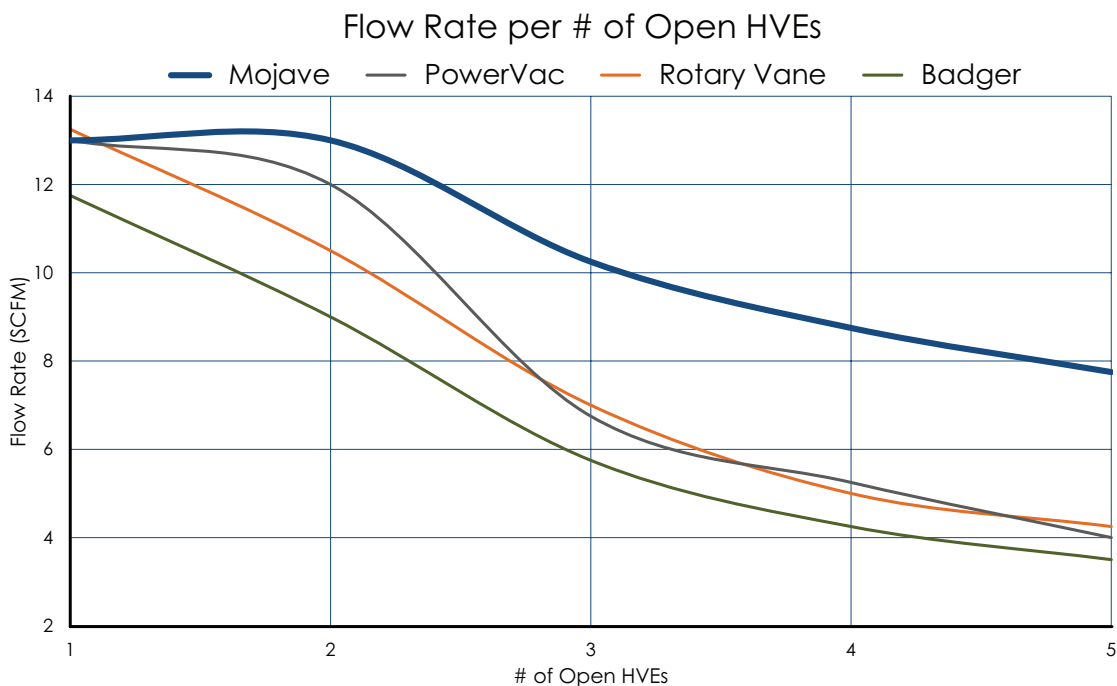


Mojave provides the highest suction power of any dry vacuum system available today. No matter how many users in your practice, Mojave delivers more flow than the competition.

## The Smart Choice

Mojave dry vacuum systems are designed to provide the highest operational performance with low operating costs and virtually no maintenance, all in a compact footprint. Whether it's time to replace your current vacuum, considering the refurbishment of your practice or even building a new facility, Mojave is your solution.

## Essential Elements For A Longer Life

- Sound-reducing impeller design.
- 3-phase motor generates higher torque and less heat.
- Extended heat-resistant motor bearings.
- All wetted metal parts are nickel plated or stainless steel.
- Aluminum high performance air-exhaust heat exchanger.
- Rugged metal chassis.

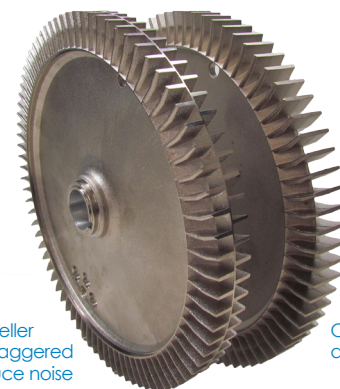

## Variable Frequency Drive

Mojave's highly efficient regenerative pumps use our market-leading Variable Frequency Drives to deliver all the power and suction you need – but only when you want it!

- Pumps speed up as users come online to maintain constant suction.
- As HVEs and SEs are shut off, the pumps slow down, reducing power consumption.
- With multiple pump systems, Mojave distributes the workload evenly.



Single Separation Tank



Dual Separation Tank

## The Mojave Difference

Mojave uses a robust 10 gallon Separation Tank designed to capture evacuate and prohibit the ingress of fluids to the pump. A single stand-alone system is perfect for smaller practices while Mojave's Continuum Tank offers uninterrupted vacuum operation for even the largest of clinics.

- Tank drains evacuate when the pumps are turned off or the internal float switch is triggered.
- Optional water connection to take advantage of the Auto Tank Washout Cycle.
- With a Continuum Tank, one tank drains as the other tank fills.
- Rugged metal chassis securely mounts all components.
- Stackable with pump chassis for a fast and easy installation.

## Mojave Master Controller

The Mojave Master Controller (MMC) automatically adjusts the speed of the pump(s) to maintain vacuum level regardless of how many users come on-line. A vacuum transducer senses when a HVE or SE is opened, increasing the frequency of the 3-phase motor to quickly provide the flow rate needed on-demand.



- Energy efficient system that prolongs the life of your pumps.
- Monitors and displays pump frequency, vacuum level and temperature.
- Set vacuum level from 8 to 10 inHg (V15 models up to 15 inHg).
- Records the run time in hours, of up to 4 connected pumps.
- Balances the vacuum load equally across all running pumps.
- Displays any error/fault codes.
- Controls the washout and diagnosis functions.



## V3 - TYPICAL PRACTICE



## 2V15 - LARGE CLINICS



# For The Smallest Of Practices, To The Largest Of Clinics

Mojave is designed with high performance, high efficiency and low operational costs. Our V3, V5 and V7 units are perfect for offices with just 3-7 simultaneous users. Choose between the economics of a single pump, or gain peace-of-mind with added redundancy from a dual-pump system. A small footprint enables easy installation in almost any facility. Big performance ensures the demands of every user are met.

The Mojave 15 series is designed for large practices, clinics and dental schools. A single pump can handle up to 15 simultaneous users and the system is fully expandable up to four parallel units supporting 60 simultaneous users.

- Highest flow performance vacuum.
- Dramatically reduces power consumption.
- Virtually eliminates water use.
- Easy installation, maintenance and expansion.

***“Mojave dry vac was an easy choice for me. It consumes very little electricity and yet provides optimum performance. Air Techniques has designed a top-rate dry vacuum system in Mojave. In addition, my dental assistant loves the innovative self-rinsing system.”***

*- Matthew Hyde, DDS*

***“When I first built my office I had 3 operatories. Later we expanded up to 10 ops. I was able to add a second pump and a continuous tank without missing an hour of practice time. The expansion was so quick our vacuum was only down for 15 minutes during lunch. It was so easy!”***

*- Richard E. Bonati, DDS*

# Specifications

Part No.	Max. Users	Horsepower	Panel Circuit Breaker	Voltage (+/- 10%)	Pump Dimensions W x H x D (in)	Tank Included	Tank Dimensions W x H x D (in)	Weights (lbs)
V3	3	1.5	20(A)	220	25 x 21 x 17	MT10	29 x 22 x 33	Pump - 145 Tank - 75
V5	5	1.5	20(A)	220	25 x 21 x 17	MT10	29 x 22 x 33	Pump - 145 Tank - 75
V7	7	2.0	30(A)	220	25 x 24 x 17	MT10	29 x 22 x 33	Pump - 145 Tank - 75
2V3	6	3.0	20(A) (x2)	220	25 x 21 x 34	MT10	29 x 22 x 33	Pump - 145 (x2) Tank - 75
2V3CT	6	3.0	20(A) (x2)	220	25 x 21 x 34	CT20	29 x 23 x 42	Pump - 145 (x2) Tank - 150
2V5	10	3.0	20(A) (x2)	220	25 x 21 x 34	MT10	29 x 22 x 33	Pump - 145 (x2) Tank - 75
2V5CT	10	3.0	20(A) (x2)	220	25 x 21 x 34	CT20	29 x 23 x 42	Pump - 145 (x2) Tank - 150
2V7	14	4.0	30(A) (x2)	220	25 x 21 x 34	CT20	29 x 23 x 42	Pump - 145 (x2) Tank - 150
3V5	15	4.5	20(A) (x3)	220	25 x 21 x 51	CT20	29 x 23 x 42	Pump - 145 (x3) Tank - 150
4V5	20	6.0	20(A) (x4)	220	50 x 21 x 34	CT20	29 x 23 x 42	Pump - 145 (x4) Tank - 150
V15	15	6.2	40(A)	220 - 3Φ	33 x 26 x 24	CT12	25 x 23 x 42	Pump - 250 Tank - 175
2V15	30	12.4	40(A) (x2)	220 - 3Φ	33 x 26 x 47	CT12	25 x 23 x 42	Pump - 250 (x2) Tank - 175
3V15	45	18.6	40(A) (x3)	220 - 3Φ	65 x 26 x 47	CT12 & CT34	50 x 23 x 42	Pump - 250 (x3) Tank - 175 (x2)
4V15	60	24.8	40(A) (x4)	220 - 3Φ	65 x 26 x 47	CT12 & CT34	50 x 23 x 42	Pump - 250 (x4) Tank - 175 (x2)
Master Controller ♦	N/A	N/A	15(A)	120* 220**	11 x 3 x 10	N/A	N/A	13

♦ One Master Controller for each system

♦ 220V Master Controller available (MT12 & CT22)

\*MMC model for V3, V5 & V7 systems


\*\*MMC model for V15 systems

## Easy Maintenance

Mojave has been designed to make maintenance easy with no filters to replace or traps to drain.

We recommend the use of Monarch CleanStream® on a daily basis to break down and eliminate the debris build-up that occurs in vacuum lines. It is the only cleaner we recommend using with Mojave Dry Vacuum Systems.



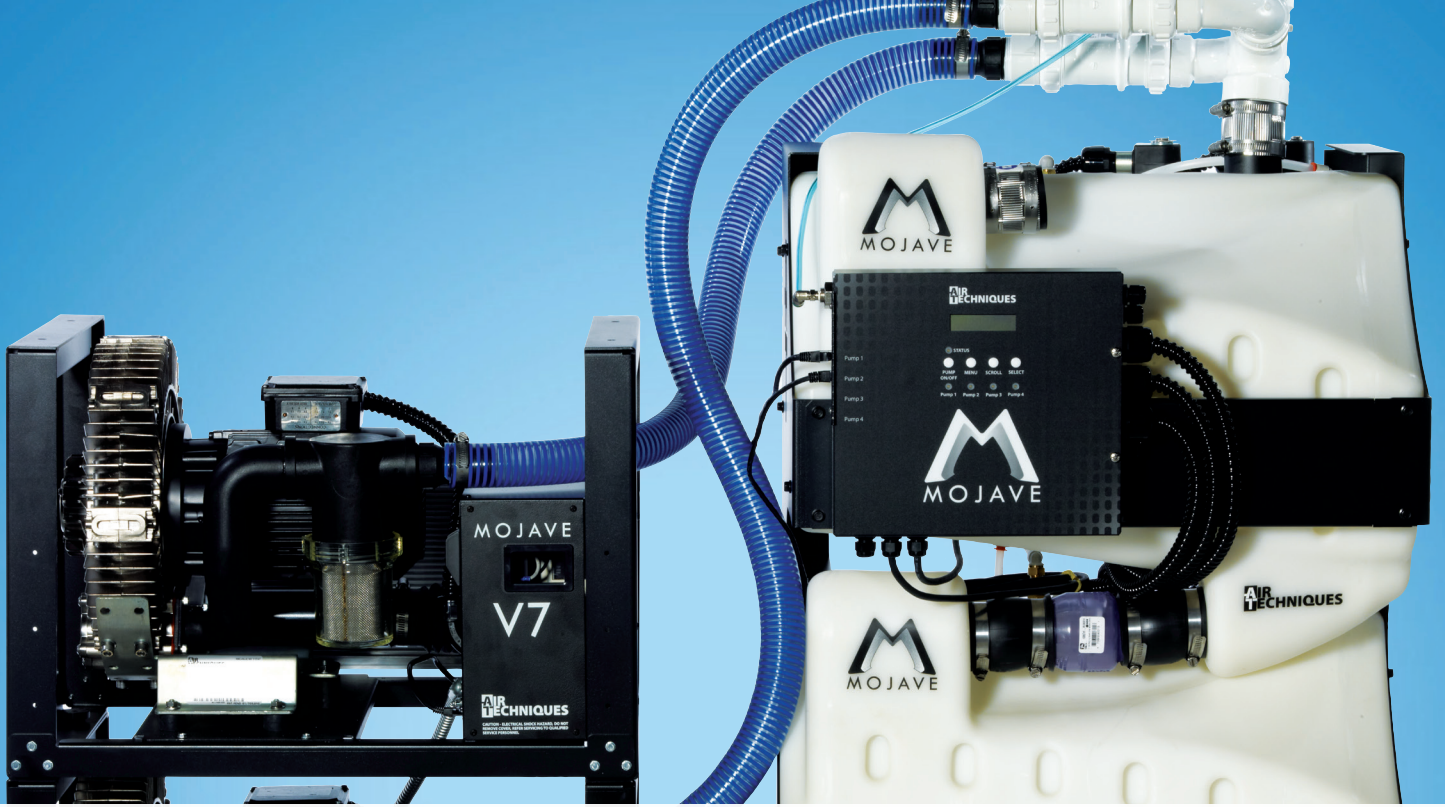

All Mojave Models Are Backed By A 5 Year Limited Warranty\*

\*Mojave motors are under warranty for 5 years

[airtechniques.com](http://airtechniques.com)

## Over 50 Years Of Proven Performance

More than 50 years ago, Air Techniques introduced its first dental air compressor. Today, Air Techniques air compressors and vacuum systems are the standard for dental utility equipment. They are used to power more dental offices in North America than all other air compressors and vacuum systems combined.

Nearly two decades ago Air Techniques introduced one of the first dry vacuum systems in dentistry designed to reduce water usage by dental offices in response to the rapidly growing costs of water. Mojave is the next generation of dry vacuum systems from Air Techniques. Mojave Dry Vacuum Systems enable dental facilities from the smallest solo practitioner to the largest of clinics experience the suction power, reliability and eco-friendly benefits only dry vacuum systems can provide.

**Air Techniques is the leader in dental air compressor and vacuum systems. Nearly 7 of 10 air compressors and 6 of 10 vacuum systems sold today are made by Air Techniques.**



CORPORATE HEADQUARTERS  
1295 Walt Whitman Road  
Melville, New York 11747- 3062  
800-247-8324



Mojave is a registered trademark of  
Air Techniques, Inc.  
© Copyright 2017  
PN 9927-775 Rev. K

**AIR**  
**TECHNIQUES** equipped for life®