

Aerosol Reduction

DuraMax[®] M800-NS

Aerosol Reduction High Speed Handpiece



Our award winning M800 direct-LED handpiece is now coupled with Safeguard Air Technology to minimize aerosol production and to reduce airborne contamination during dental procedures.

The DuraMax M800 has the same powerful 23W of cutting power and brilliantly bright illumination without fiber-optic system needed, but with increased visibility during operation and enhances patient and office safety.

Make stress-free treatment a reality.




How can this handpiece reduce aerosol generation?



STREAM WATER

This handpiece is equipped with single stream water to cool down the bur instead of traditional spray water.

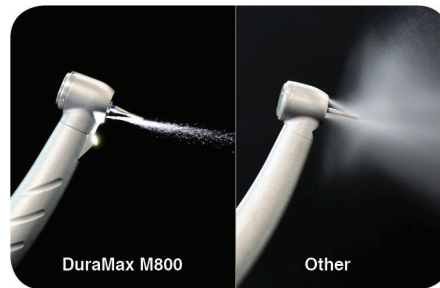

SAFEGUARD AIR[®]

There are three specialized arrayed ports around the bur to push aerosol downward.

Key Specification

Illumination Source	Built-in Generator and Direct LED
Colour temperature(K)	5,500
Brightness(LUX)	25,500
Head Size(Width x Height)	10.8mm x 13.0mm
Drive Air Pressure(PSI)	40-60 PSI
Power(WATT)	23
Speed(RPM)	410,000
Rotary Instruments	FG burs according to EN ISO 1797-1:1995
Ø / max. length(MM)	1.6 / 25
Chuck Type	Push Button
Water Port	Single Stream Water

Order Information

Beyes® PD Connection



EZ-Swivel® PD Coupler



Code **HP3056**
 Model **M800-M/PD/NS**



Code **HP6008**
 Model **PD M4**



Beyes®
 The Next Standard

Beyes Dental Canada Inc.
 595 Middlefield Road
 Toronto, Ontario, M1V 3S2
 Canada
 Tel: 1-416-281-1888
 Email: info@beyes.ca
 Website: www.beyes.ca

Your Dealer for Beyes Products

Scan the QR code for more information



Printed in Canada
 File #: ENB029
 Rev1/23.06.20

© Beyes Canada Inc 2020. All Rights Reserved.
 Specifications are subject to change without notice.
 All brands are holders of their respective trademarks.
 All Beyes products are manufactured in strict compliance with Health Canada regulations ensuring quality.
 Manufacturing facility is ISO 13485 certified under Canadian Medical Devices Conformity Assessment System (CMDCAS).