

# Maxso<sup>®</sup>

## High Speed Handpiece

The Maxso series high speed handpiece is designed to meet expectations in modern dental procedures. High quality components and comfortable grip with Micro-tex finish provide superior handling that you desire. Furthermore, the well aimed triple water spray will provide excellent coolant effects. Maxso series comes with variety of connections; as a result, you simply attach the Maxso to your existing system.



### Triple Water Spray

Effectively cools the entire surface of the bur that avoids overheating and provides better patient comfort.





### Micro-Tex Coating

The latest surface coating technology creates a micro-textured and scratch proof to enhance handling and durability.



### Compatibility



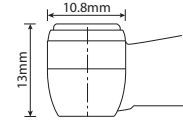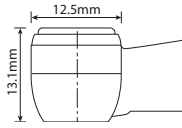


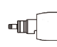


Compatible with a great variety of coupling systems and adapts great to every existing treatment units.



## Technical Data

### Standard Head

### Mini Head

	Standard Head	Mini Head
<b>Head Size(Width x Height)</b>	12.5mm x 13.1mm	10.8 x 13.0mm
<b>Drive Air Pressure(PSI)</b>	40-45 PSI	40-45 PSI
<b>Power(WATT)</b>	21	19
<b>Speed(RPM)</b>	350,000	400,000
<b>Rotary Instruments</b>	FG burs according to EN ISO 1797-1:1995	FG burs according to EN ISO 1797-1:1995
<b>Ø / max. length(MM)</b>	1.6 / 25	1.6 / 25
<b>Chuck Type</b>	Push Button	Push Button
<b>Water Spray</b>	Triple	Triple
<b>Midwest 4 Hole Connection</b>		
 	CODE <b>HP2008</b> MODEL <b>M200-S/M4</b>	CODE <b>HP2009</b> MODEL <b>M200-M/M4</b>
<b>Borden 2 Hole Connection</b>		
 	CODE <b>HP2012</b> MODEL <b>M200-S/B2</b>	CODE <b>HP2013</b> MODEL <b>M200-M/B2</b>
<b>Beyes® PD Connection</b>		
 	CODE <b>HP6022</b> MODEL <b>M200-S/PD</b>	CODE <b>HP6012</b> MODEL <b>M200-M/PD</b>
<b>NSK® Connection</b>		
 	CODE <b>HP2036</b> MODEL <b>M200-S/N</b>	CODE <b>HP2037</b> MODEL <b>M200-M/N</b>

**Beyes Dental Canada Inc.**  
 30 Casebridge Court  
 Toronto, Ontario, M1B 3M5  
 Canada  
 Tel: 1-416-281-1888  
 Fax: 1-855-720-1228  
 www.beyes.ca  
 info@beyes.ca

**Beyes EC REP**  
 26 York Street  
 London, W1U 6PZ  
 United Kingdom  
 Tel: 44-20-3514-0949  
 Fax: 1-855-720-1228  
 www.beyes.ca  
 info@beyes.ca



C108EN.REV0E1515  
 PRINTED IN CANADA

\*All brands are holders of their respective trademarks.