

DENTAPEN

by Septodont

the **PERFECT MATCH**
for your practice

NEW!



DENTAPEN

Beyond injection: Dentapen, the new generation of powered injector for Dental Anesthesia

Matches your patients' preferences: constant flow to reduce pain, non-threatening look

Matches your needs: comfort and reproducibility of your injections all-day-long

Matches your habits: no training, both syringe and pen-like handling, compatible with any needle

MANAGING
PAIN FOR
YOUR
PRACTICE



CCLAD TECHNOLOGY*

TESTED ON DENTISTS

● A study conducted on dentists demonstrated a marked reduction in pain perception for injections using a CCLAD system**


● **96%**  of dentists tested received a palatal injection, and preferred the electronic syringe over the manual one

● **Pain was reduced** when compared to a standard manual syringe



**CCLAD (Computer-Controlled Local Anesthesia Device)*

***Hochman MN, Chiarello D, Hochman CB, Lopatkin R, Pergola S. Computerized Local Anesthesia Delivery vs. Traditional Syringe Technique. NY State Dent J. 1997;63:24-9.*

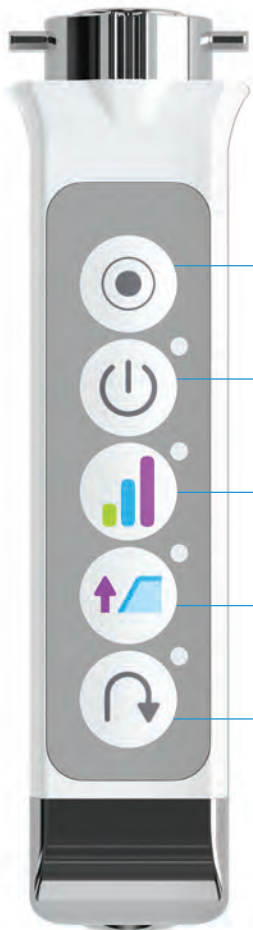

NO TRAINING NEEDED

5 INTUITIVE BUTTONS

- **Activation**
automatic self-aspiration
- **On / Off**
with battery level indicator
- **Speed Selection**
fast: 1ml/30s | medium: 1ml/60s | slow: 1ml/90s
- **Mode Selection**
 - intraligamentary mode: the slowest speed with maximum force
 - ramp-up mode: anesthesia flow increases gradually and constantly
- **Rewind Retraction**
return to the initial position and insert the next cartridge

EFFORTLESS INJECTION

That matches your preferences



Ergonomic
and light (40g/1.4oz)



Compatible
with the needle of your choice



Standard syringe grip
or Nerve block finger grip



SYRINGE-LIKE HOLDING



PEN-LIKE HOLDING

PRACTICAL

Focus your attention on the injection



No cord
minimizes pressure on hand
and needle



No console
space saving



No foot activation
focus attention on the injection

FOR YOU



Shown actual size.

FOR YOUR PATIENTS

A MORE COMFORTABLE INJECTION

Improves your patients' satisfaction with less fear of injection



Constant flow
2 modes



3 speeds



Small and modern
less fear of injection

THE WORLD LEADER IN PAIN MANAGEMENT

Septodont has always been at the forefront of dentistry, providing you with a wide range of products and services to help you administer a painfree experience.

Septodont continues to innovate with the Dentapen - the perfect match for your practice.



DENTAPEN KIT INCLUDES

3x
Batteries

1x
Dentapen Body

5x
Finger Grips

7x
Cartridge Holders

25x
Protective Sleeves

Accessories

#01N6010: Refill Finger Grips Syringe Style (5 per box)

#01N6020: Refill Finger Grips Pen Style (5 per box)

#01N6030: Refill Cartridge Holder (5 per box)

#01N6040: Refill Battery (10 per box)

#01N6050: Refill Protective Sleeves (250 per box)

All accessories are available separately.
Finger grips and cartridge holders must be sterilized before use.



#01N6000: Dentapen Kit

MANAGING PAIN FOR
YOUR PRACTICE
www.septodontusa.com

